



Liquid Balance

CS5 supports you to collect the infusion volume automatically. By management of other input and output flow volume, allows you to integrate the data into one liquid balance graph. Help nurses and doctors to easily master patient's condition, and make further treatment solution.



History Record Management

CS5 allows you to monitor alarm information system of all wards. With the powerful storage ability, CS5 supports for at least 10 days/2000 pieces record per patient, 20,000 patients in total. With the classification of event type, allows you to easy track for event review, treatment review, as well as drug trend review.

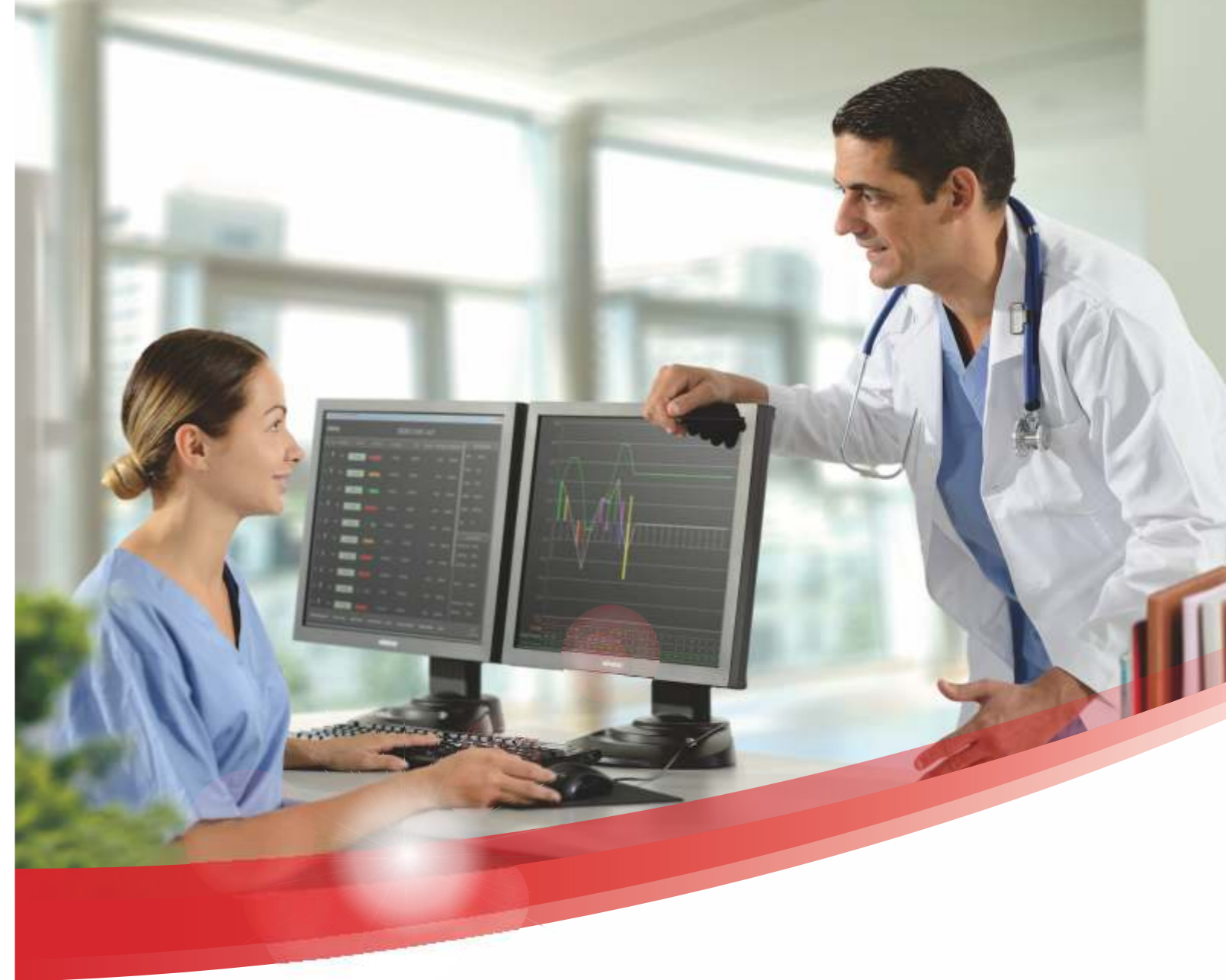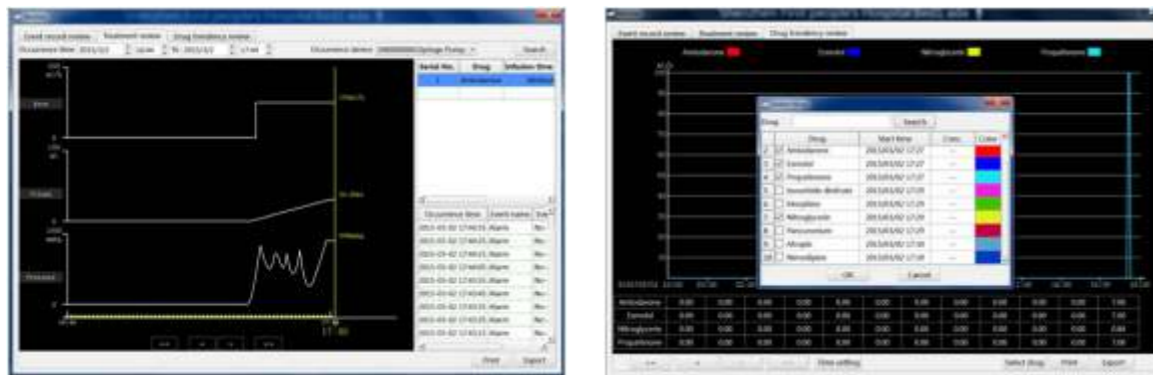

BeneFusion CS5

Central Infusion Supervision System

Shenzhen Shenke Medical Instrument Technical Development Co., Ltd.
 Floor 2, 3, 4, Building 1, Section 1, Baiwangxin Hi-Tech Industrial Park, Songbai Road, Xili, Nanshan, Shenzhen 518108, P.R. China.
 E-mail: skmedica@skmedica.com
 Website: www.skmedica.com


mindray is a trademark of Shenzhen Mindray Bio-Medical Electronics Co.,Ltd.
 ©2013 Shenzhen Mindray Bio-Medical Electronics Co.,Ltd. All rights reserved. Specifications subject to changes without prior notice.
 P/N:ENG-BeneFusion DS5 -210X285-20170619 046-007839-00(1.0)

mindray

mindray
 healthcare within reach


BeneFusion CS5 Central Infusion Supervision System

BeneFusion CS5 central infusion supervision system allows you to central monitor the infusion status and real-time information of all patients and medical pumps. The powerful operation system allows for connecting with at least 64 units of BeneFusion DS5 infusion supervision system, and over 512 units of infusion/syringe pumps.

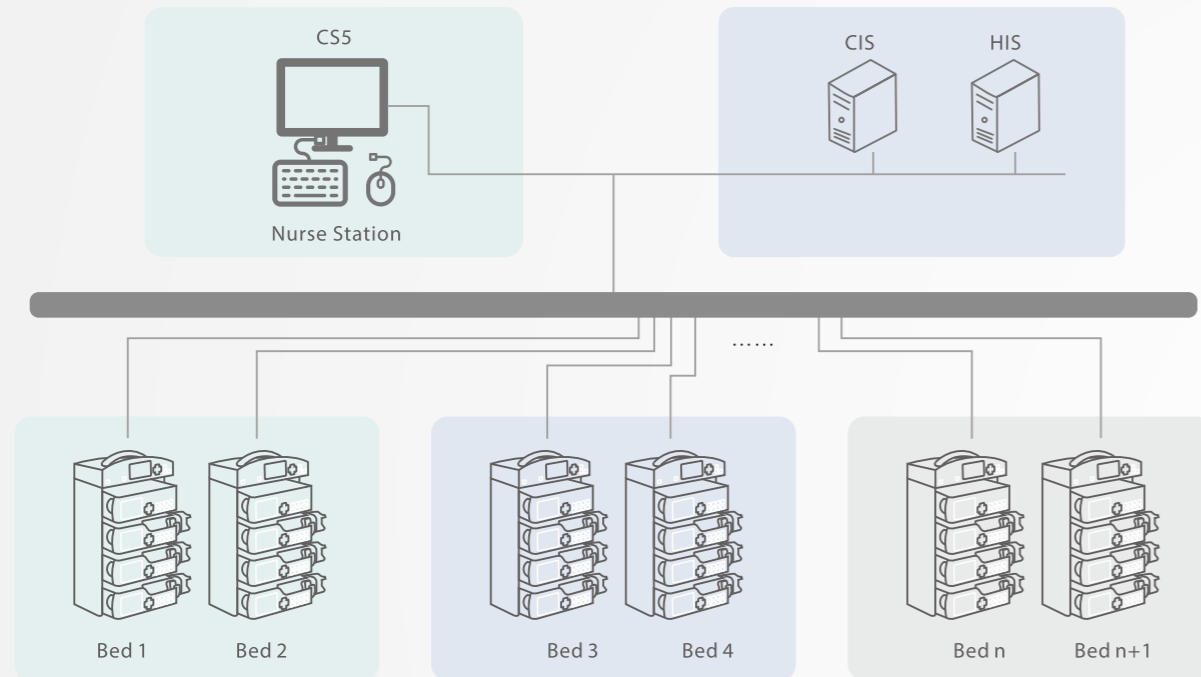

Patient and Infusion Information Management

Through BeneFusion CS5 central infusion supervision system, allows you to draw and display the actual location of any beds in ward, helping you to deal with clinical problems quickly. Admit patient and discharge patient, real-time display, check and edit patient information manually, or through barcode scanner, support at least 20,000 patients infusion information, alarm event and treatment record view.

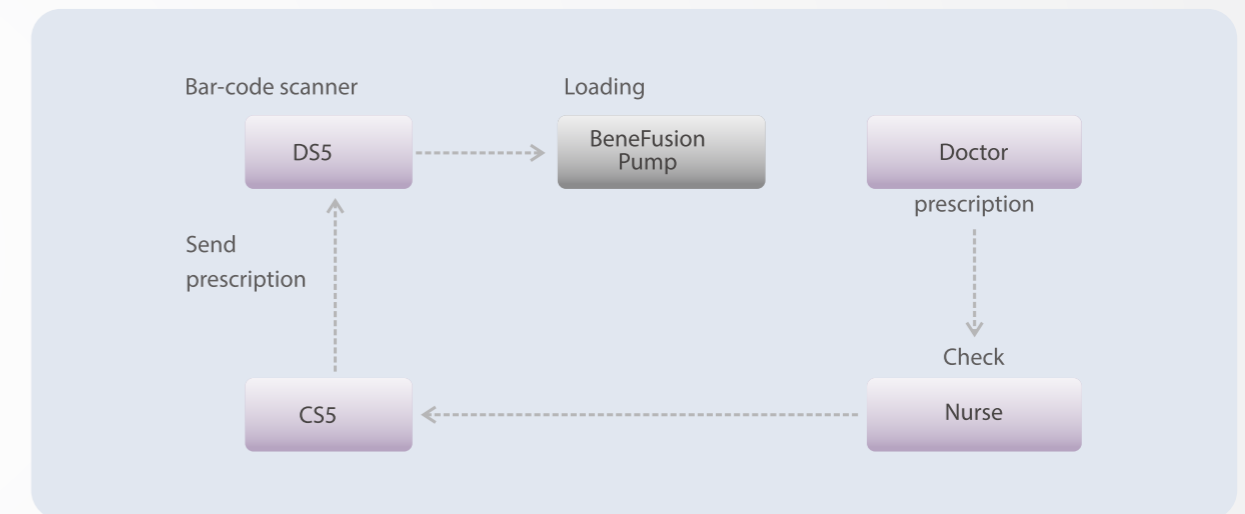

CS5 could automatically download the prescription from CIS/HIS, then send it to DS5 directly. Nurses do not need to edit the drug information any more, just press the start button, the pump will start the infusion quickly as the loaded prescription.

Drug Library Management

CS5 allows you to edit/import/export drug library, and easy updating for BeneFusion DS5 or infusion/syringe pump directly. The general drug library version supports for at least 2000 drugs.

