



**BeneHeart D1**  
Defibrillator





## BeneHeart D1

Defibrillator



## Powerful support for professional responders

BeneHeart D1 is a powerful automated external defibrillator (AED). Supporting you throughout the emergency response process, BeneHeart D1 is lightweight, compact, durable and very easy to use while offering all standard AED functions.

Based on biphasic technology, the delivered energy can be escalated to 300J and 360J. BeneHeart D1 allows you to improve the defibrillation success ratio by increasing the dose of subsequent shocks.

Acknowledging the stressful and time-sensitive working environment, D1's professional grade configuration has been designed especially for **BLS** and **ALS** teams to treat cardiac arrest patients and provide continuous ECG monitoring. The intelligent functions and intuitive user interface make D1 an **efficient** and **flexible** lifesaving aid for professional responders.









## Efficient: Helping BLS to respond faster

- **Large 7-inch TFT display** supports clear visual guide with color animation.
- **Pre-connected pads** save time.
- **Easily switch between pediatric and adult** patients with the press of one button.
- **Auto-identify patient type** through pad connector so that a lower energy dose is delivered to pediatric patients.
- **Auto-adjustable volume and screen brightness** allow you to comfortably use D1 in a variety of environments.

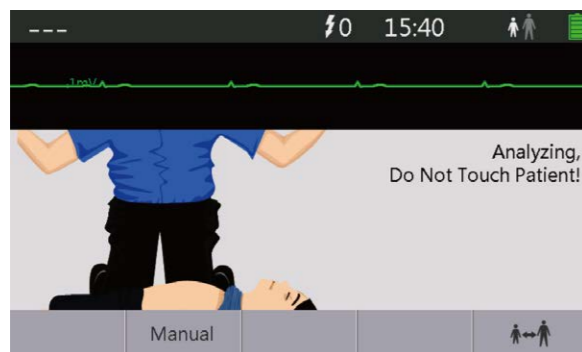

## Flexible: Helping ALS to provide more professional therapies

- **Real manual mode with energy level selection** provides the complete 1-2-3 step process for you to deliver high-quality shocks to the patient.
- **Independent 3-lead ECG monitoring** through separate socket allows you to view more ECG leads or arrhythmias while the pads are still being used.
- **Synchronized cardioversion** is also available in the small and compact D1.

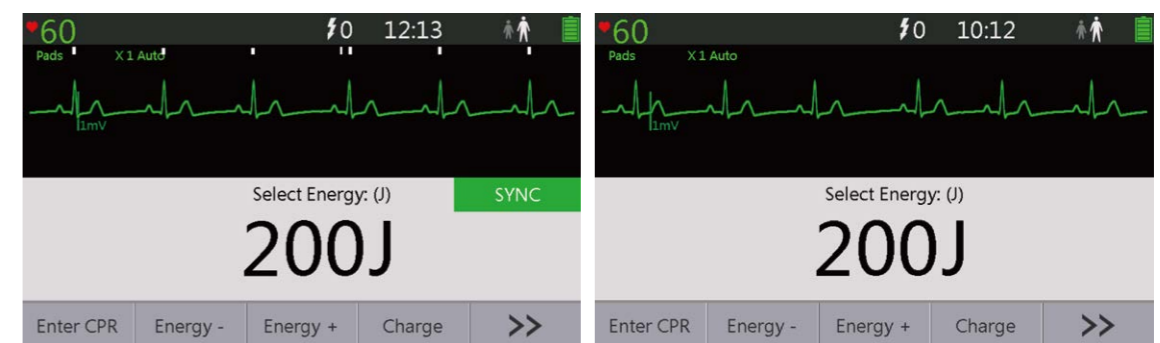



# BeneHeart D1

## Defibrillator

### Battery and External Charger



BeneHeart D1 with one Li-ion battery supports more than 300 shocks at 200J or 200 shocks at 360J as well as 12 hours of monitoring. The external charger keeps a spare battery charged and ready for use.

\* The Li-ion battery is compatible with BeneHeart D3's.

### LED Indicator



Automated self-tests are performed periodically and results are shown through audiovisual alerts. Also, a visual indicator in green confirms that D1 is ready for use.

### Carrying Case



The carrying case offers space for D1 together with a spare battery, ECG wires and cable as well as a small kit containing scissors and a razor.

### USB Port



You can export data from D1 through a USB drive and review it on a PC using Mindray software.

## CPR Solution improving compression quality intuitively



### According to the 2015 AHA/ERC guideline

The CPR Sensor is easy to use with its ergonomic design and graphic mark of the pressure spot. It delivers instant audiovisual feedback so that compressions are of proper depth and rate. It indicates complete chest recoil, hands-off time to help improve the patient's chance to survive.

CCI (CPR compression index): display dynamic change of CPR compression quality, including information such as compression rate, compression depth, compression recoil, compression fraction

Compression waveform area

Compression rate area & Compression depth area

